


Full Face Shield / Eye Protection Reuse During Contingency Level

Effective: 3/26/2020

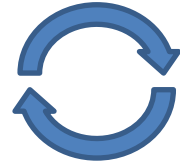

Follow these guidelines to re-use full faces shields, goggles, and disposable faces shields when indicated

Discard your Full-Face Shield /Eye Protection if any of these apply

- The face shield or foam is visibly soiled, torn, damaged, or can no longer fasten securely.
- The face shield is difficult to see through



Otherwise, save your face shield; to be reused by the same healthcare provider. It may be used for encounters with different patients.

Doffing a Used Face shield

1. Perform hand hygiene
2. Put on gloves
3. Remove face shield by touching the elastic band only
4. Disinfect the outside of the face shield/eye protection with EPA approved disinfectant wearing clean gloves
5. Wipe the outside of the face shield/eye protection with water or alcohol to remove residue
6. Air dry or dry using clean absorbent towels
7. Drop in large paper bag
8. Remove gloves, perform hand hygiene
9. Label the bag with your name
10. Re use face shield/eye protection for the duration of shift.
11. Fold Paper bag & store in designated unit location for your shift. Discard face shield/ eye protection at end of shift.



Donning a Used face shield

1. Grab paper bag & open
2. Perform hand hygiene before removing the face shield from stored location.
3. Verify that the face shield is not visibly soiled, torn, nor damaged
4. Don the face shield.
5. Once the face shield is donned, perform hand hygiene